

The logo for 'enova' features a stylized green and blue 'e' above the word 'enova' in a blue sans-serif font.

enova®



**BROKEN INSTRUMENT REMOVAL KIT  
& MICRO TWEEZER**


The 'zumax' logo consists of the word 'zumax' in white lowercase letters inside a blue rectangular box.

zumax

# Broken Instrument Removal Kit (PN FT0017-25)



## Kit Sizing Options

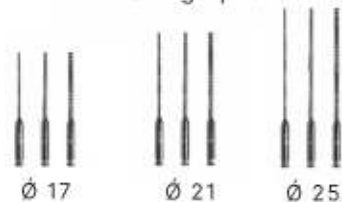

## Extractor Tip Options



## Advantages:

- Limits amount of dentin removed from the root canal wall.
- Reduces risk of the ultrasonic tip fracture in a traditional technique & secondary fracture of the separated nickel titanium instrument when using the ultrasonic technique.
- Microscopic tweezers are favorable and efficient for safe operations.
- The operation process is quick and easy to learn.

## Contents:

The broken instrument removal kit contains 3 micro-trepans, 1 extractor handle, 3 extractors, 1 foreign forceps, and 1 instrument box.

Extractors are mounted on the extractor handle for use.




# Micro Dental Tweezer

(PN DT0001)

## Instructions for Dental Tweezers:

1. When any foreign matter is discovered during procedure, insert tweezer tip into canal to remove matter.
2. After foreign matter is removed, use the cleaning adapter to clean interior of tweezer.
3. Autoclave entire tweezer and tip after each use.

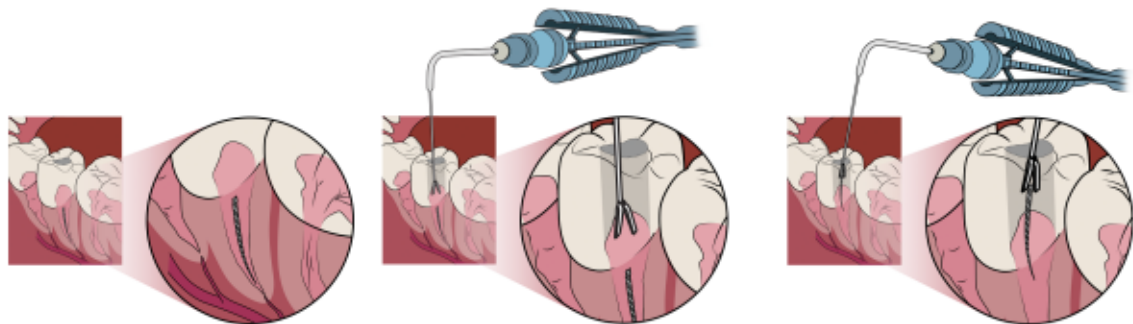

## Advantages:

- 360 degree adjustment.
- Head shape is designed for easy clamping.
- Accurately clamps foreign material with firm grip, high stability, and reliability.

## Contents:

Each package includes 1 tweezer head, 1 handle, and 1 cleaning adapter.  
Specifications: 0.9mm x 20mm



Micro-Claw  
(0.9mm x 20mm)

To learn more about how Enova can help you see better, so you can perform better, visit [www.goenova.com](http://www.goenova.com), call 651.236.8857, or email [sales@goenova.com](mailto:sales@goenova.com).

**World Headquarters**

Enova Illumination  
1717 Broadway Street NE  
Minneapolis, MN 55413



© 2021 Enova Illumination, Inc. All trademarks are the property of Enova Illumination, Inc. unless otherwise noted. F20414

Zumax™ is a Trademark of Zumax Medical Co., Ltd.